



Condensation particle counter for number concentration up to  $10^4$  particles/cm<sup>3</sup> (single count mode) and  $10^7$  particles/cm<sup>3</sup> (photometric mode)

## Description

The Palas® **UF-CPC 50 condensation particle counter** is optimized for low concentrations (e.g. downstream from a filter medium or in operating rooms).

The UF-CPC 50 measures the total particle concentration of ultrafine and nanoparticles suspended in air or other carrier gases. These particles are enlarged by a condensation process in order to enable the precise determination of their number with an optical light scattering detector. The working fluid, e.g. butanol or water, is vaporized as a condensation agent. The nanoparticles to be measured are directed through the vapor atmosphere, during which the vapor condenses on the nanoparticles in a cooling zone. The condensation process is influenced by the nanoparticles themselves, as well as the working fluid, the operating temperatures, and the volume flow.

In addition to their number, the UF-CPC 50 also measures the size of the droplets and provides the user with additional information concerning the condensation process. This novel and patented (German patent no.: DE<tcxspan> 10 2005 001 992 </tcxspan>A1) working fluid delivery system provides greater flexibility. The user is no longer limited to butanol or water but may also opt for more environmentally friendly or better suited working liquids for special applications, in addition to the butanol and water.

In addition, the modularity of the components enables the user to perform most maintenance tasks (e.g. cleaning, pump replacement).

In research mode, the user is able to easily adjust a variety of parameters, e.g. temperature settings on the saturator, using the 7" touch screen.

For process control applications, the UF-CPC supports a standardized interface with various protocol options, e.g. Modbus, and provides features, such as remote access and data storage on the Internet or an internal network.

Figure 1 presents the principle of operation of the UF-CPC [1]. The aerosol with nanoparticles enters the UF-CPC at the bottom and reaches the heated evaporation chamber – the saturator first.

Within the saturator, the working fluid is moved helically around the flow area of the aerosol, resulting in a more homogeneous contact area compared to designs where only one or two walls of the saturator are lined with a porous material that is soaked with the working fluid

In addition, the working fluid is circulated continuously from the reservoir to the constantly heated helical and U-shaped channel and back to the reservoir at a flow rate that can be adjusted to accommodate different working fluids.

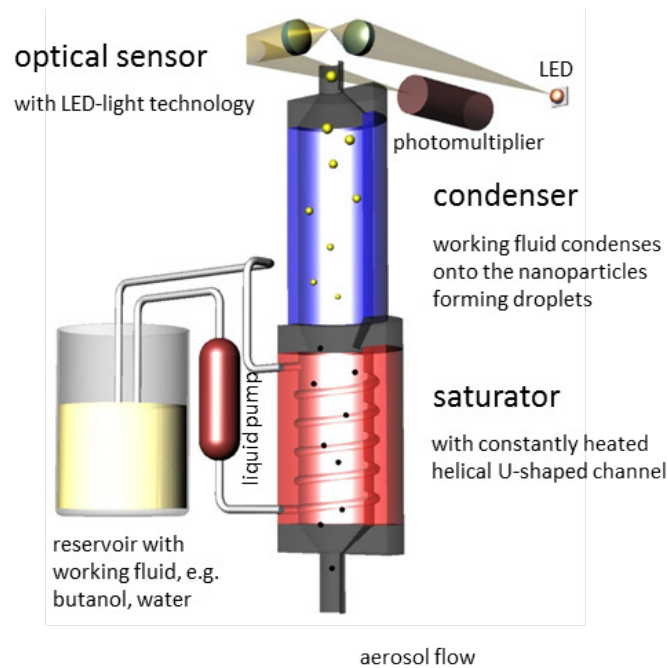


Fig. 1: Principle of operation of the condensation particle counter (UF-CPC)

Downstream from the evaporation chamber, the aerosol and the saturated carrier gas enter a cooled area, the condenser, in which the working fluid condenses onto the nanoparticles, forming droplets in the  $m$  size range.

Downstream from the condenser, the droplets enter the optical sensor. Here, the size of the droplets is analyzed and their concentration is measured by counting them. Unlike other CPCs, the sensor on the UF-CPC uses a patented technology to count particles at concentrations of up to  $10^7$  P/cm<sup>3</sup> without diluting the aerosol.

### Software

Based on continuous customer feedback, the user interface and software have been designed for intuitive operation and real-time control of measurement data and parameters (Fig. 2) In addition, the software provides data management with the integrated data logger, sophisticated export capabilities, and network support. The measured data are able to be displayed and evaluated with many available options. If a specific display is desired, we can implement it for you.

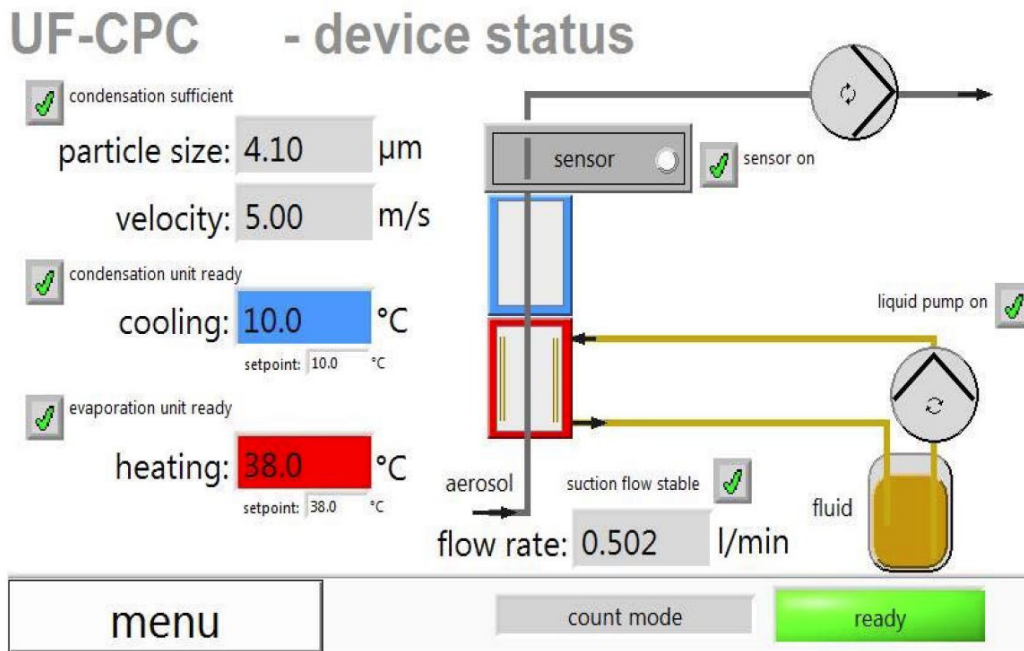


Fig. 2: Touch screen user interface showing the device status

### References

[1] Patent: DE 10 2005 001 992 A1

## Benefits

- This novel and patented working fluid delivery system provides the user with the option to quickly switch between butanol and isopropanol, for example, or even water.
- The UF-CPC is able to count up to 2,000,000 particles/cm<sup>3</sup> in single count mode, depending on the sensor in use (able to be switched out by the user).
- Integrated computer with 7" touch screen
- Intuitive user interface with sophisticated software for data analysis
- Integrated data logger
- Unlimited network compatibility that supports remote control and data storage on the Internet
- Integrated interface for process control applications

## Datasheet

Parameter	Description
Measurement range (size)	4 - 10,000 nm
Measurement range (number C <sub>N</sub> )	10 <sup>4</sup> particles/cm <sup>3</sup> (single count mode), 10 <sup>4</sup> - 10 <sup>7</sup> particles/cm <sup>3</sup> (nephelometric mode)
Volume flow	0.9 l/min (butanol); 0.3 - 1 l/min (adjustable for research applications)
Data acquisition	Digital, 20 MHz processor, 256 raw data channels
Light source	LED
User interface	Touch screen, 800 • 480 pixel, 7"
Dimensions	29 • 24 • 35 cm (H • W • D)
Data logger storage	4 GB
Software	PDAnalyze
Accuracy	5% (single count mode); 10% (nephelometric mode)
Response time	t <sub>90</sub> = 2.8 s, t <sub>90-10</sub> = 2.0 s
Operation liquid	Butanol (standard); for research applications, further liquids as isopropanol or water can be used and exchanged easily
Installation conditions	+10 - +30 °C (others on demand)

## Applications

- Aerosol research
- Test of filters and air cleaners
- Environmental measurements
- Studies on working place pollution and working place safety
- Studies on inhalation and health effects
- Process surveillance
- Studies on emission of printers

**Palas GmbH**  
Partikel- und Lasermesstechnik  
Greschbachstrasse 3 b  
**76229 Karlsruhe**  
Germany

**Managing Partner:**  
Dr.-Ing. Maximilian Weiß, Udo Fuchslocher  
**Commercial Register:**  
register court: Mannheim  
company registration number: HRB 103813  
USt-Id: DE143585902



**Contact:** E-Mail: [mail@palas.de](mailto:mail@palas.de) Internet: [www.palas.de](http://www.palas.de) Tel: +49 (0)721 96213-0 Fax: +49 (0)721 96213-33